

Invacare® | **Leo™**



**Safe, solid, secure: the *Invacare Leo***

The *Invacare Leo* is a 4-wheel scooter designed for all those that value their independence and wish to get out and about unaided. Safety is a key feature of the **Leo**, but this does not detract from its stylish and sporty looks. **Leo** offers people the freedom and confidence to enjoy their essential daily outings and leisure excursions. With high standards of design and engineering safety, the **Leo** undoubtedly sets a new standard.



**Yes, you can.®**



**For a safe, powerful drive**

The *Invacare Leo* offers a smooth, safe drive in a scooter that is easy to control and manoeuvre. The **Leo** intelligent lighting package includes indicators that automatically switch off and a rear brake light for safe driving. The lighting system is designed so that if a bulb should blow, all other lights still continue to function. In addition to a visual warning for the battery level display, drivers are also alerted with an audible warning (beep) when the battery needs recharging.

**Compact**

The **Leo** is a compact mini-scooter that can be easily dismantled for transportation during longer trips or weekend breaks. With a smooth and manoeuvrable drive, getting round tight corners in the local shopping centre, museum or leisure complex becomes refreshingly easy. Access to the battery and electronics for maintenance or servicing is much easier with the **Leo** as only one part needs to be removed.

**Reliable**

Customers can trust the **Leo** to take them wherever they need to go. From a trip to the shops or a ride to the park, the **Leo** is rendered sturdy and reliable thanks to its exceptionally high quality electronic engineering and power. These elements are protected from corrosion and water,

(splash protection) which safeguard the vital workings of the scooter against dirt and water. A two-step disengaging lever protects the scooter from free-wheeling if inadvertently knocked out of drive mode.

**Features and Accessories**



**Finite tiller movement**  
Customers can easily adjust the tiller to suit their needs.



**Brake light**  
Rear brake light warns other people that the scooter is slowing down, thus preventing collisions.



**Comfort seating**  
Swivel function and sliding rails ensure an optimal seating position.



**Two-step disengaging lever**  
This unique concept prevents the scooter from accidentally freewheeling.

**Features and Accessories**



**Built in anti-splash guards**  
Protects electronics and the transaxle from dirt and water.



**Disassembly into five parts**  
No tools are required to dismantle *Invacare Leo*



**Crutch & cane holder**  
Keeps your crutch or cane stable during the drive.



**Lockable front box**  
Secure box for personal valuables.



**Seatbelt**  
Keeps a person well-positioned in the seat.



**Large mirror**  
Can be mounted left or right.



**Heavy-duty storage cover**  
Keeps the scooter clean and dry.



**Charging off-board**  
Simply remove batteries to charge with the off-board charger kit.



**Seat suspension module**  
Adds extra comfort during driving.


## Technical data

Invacare <b>Leo</b>	470 mm	410 mm	460 - 560 mm	475 mm	590 mm	1220 mm

Invacare <b>Leo</b>	59 kg (excl. batteries) 83 kg (incl. batteries)	136 kg	2 x 36 Ah	2620 mm	60 mm	70 mm

Invacare <b>Leo</b>	10°	240 W	36 km	6 km/h

## Colour of shrouds

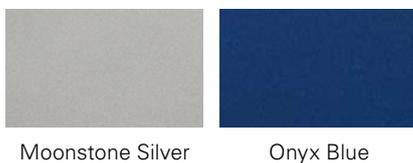

## Upholstery



Invacare Limited  
 Pencoed Technology Park - Pencoed -  
 Bridgend CF35 5AQ - United Kingdom  
 Tel: +44 1 656 776 222  
 Fax: +44 1656 776 220  
 E-mail: uk@invacare.com  
 Sales Order E-mail:  
 ordersuk@invacare.com  
 www.invacare.co.uk

Invacare Ireland Ltd  
 Unit 5 - Seatown Business Campus  
 Seatown Road - SWORDS - County  
 Dublin - Ireland  
 Tel. +353 1 8107084  
 Fax +353 1 8107085  
 Email: Ireland@invacare.com  
 www.invacare.ie

© 2011 Invacare International Sàrl  
 All rights reserved.

All information quoted is believed to be  
 correct at time of print.  
 Invacare reserves the right to alter  
 product specifications without prior  
 consultation.

Leo - UK - 08/2011 - 1503875

